

## Certificate of Performance

### TON 1

Product Code KR103

Date 02.06.2020


Expiry Date 29.05.2021

Reagent	Lot Number
TON 1	<b>4413</b>
TON 2	4414
TON 3	4415
Test Method	WC SOP 23
Test Instrument	Konelab 20
Serial Number	P4919028

Std. Value	Lot #	Abs.
0.00		0.004
2.50	TON 100 03 20	0.177
5.00	TON 100 03 20	0.342
7.69	TON 100 03 20	0.508
10.00	TON 100 03 20	0.627

\*serial dilution

### Calibration Curve

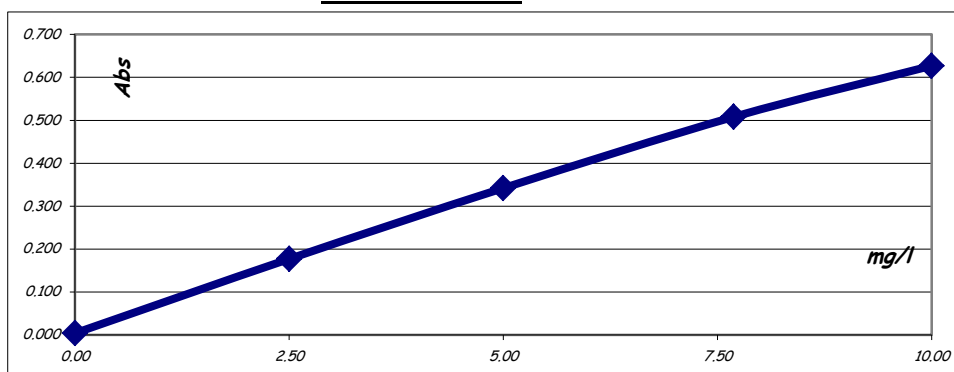

Test	Result	Specification	Pass / Fail
Intercept	0.004	-0.005 to 0.01	Pass
Sensitivity	Abs/ppm		
0 - 2.50 mg/l	0.0692	0.0550 to 0.0800	Pass
2.50 - 5.00 mg/l	0.0660	0.0600 to 0.0800	Pass
5.00 - 7.69 mg/l	0.0617	0.0550 to 0.0800	Pass
7.69 - 10.0 mg/l	0.0515	0.0500 to 0.0700	Pass

Notes

- Store between 2 & 8 °C
- \*Stable on Board TON 1 - 1 day
- \* TON 2 & TON 3 - 5 days
- Keep Tightly Sealed when not in use.
- Do not return unused portions to container
- Barcode ID- TON 1 103; TON 2 104; TON 3 105
- \* Based on results obtained in our Research Facility

*Alan Coonerty*  
A.Coonerty

An ISO 9001:2008 Registered Company

Annacotty Business Park, Annacotty, Limerick, Ireland.

Telephone 061 358190

Facsimile 061 358191

### Hazard Identification

Labelling according to Regulation (EC) No 1272/2008


 Signal Word **Warning**
**H315** Causes skin irritation

**H319** Causes serious eye irritation

 see full MSDS on [www.serosep.com](http://www.serosep.com)

Precautionary statements

P280;P264;P321

P305+P351+P338

P302+P352; P332+P313

P362+364; P337+P313