

Certificate of Performance

Ammonia 2 Reagent

Product Code **KR 109**

Lot Number **4421**

Date 10.06.2020


Expiry Date 04.06.2021

Reagent	Lot Number
AMM 1	T039
Ammonia 2	4421
Test Method	WC SOP 22
Test Instrument	Konelab 20
Serial Number	P4919028

Std. Value	Lot #	Abs.
0.00		0.019
0.25	AMM 10 06 20	0.215
0.50	AMM 10 06 20	0.415
0.77	AMM 10 06 20	0.645
1.00	AMM 10 06 20	0.833

*serial dilution

Calibration Curve

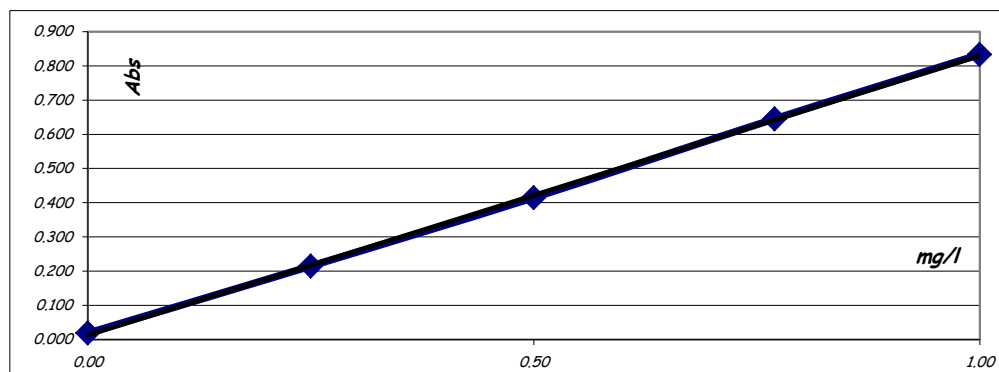

Test	Result	Specification	Pass / Fail
Intercept	0.0138	0.000 to 0.025	Pass
Sensitivity Abs/ppm	0.8167	0.750 to 0.890	Pass
Regression	0.9998	0.990 to 1.0001	Pass

Notes
 Store between 2 & 8 °c
 *Stable on Board 5 days
 Barcode ID = 109
 *Based on results obtained in our Research Facility

Alan Coonerty
 A. Coonerty

An ISO 9001:2008 Registered Company
 Annacotty Business Park, Annacotty, Limerick, Ireland.
 Telephone 061 358190 Facsimile 061 358191

Hazard Identification

Labelling according to Regulation (EC) No 1272/2008



Signal Word **Danger**
H314 Causes severe skin burns and eye damage
H318 Causes severe eye damage
H412 Harmful to aquatic life with long lasting effects

Precautionary statements
 P304+P340
 P260; P280; P310 P321; P405; P501
 P303+ P361+ P353
 P305+P351+P338
 see full MSDS on www.serosep.com