

Certificate of Performance

Fluoride 3

Product Code KR133

Date 25.06.2020


Expiry Date 24.06.2021

| Reagent | Lot Number |
|-----------------|--------------|
| Fluoride 1 | 4428 |
| Fluoride 2 | 4429 |
| Fluoride 3 | 4430 |
| Test Method | WC SOP 28 |
| Test Instrument | Gallery Plus |
| Serial Number | 862100110267 |

| Std. Value | Lot # | Abs. |
|------------|-----------------|-------|
| 0.00 | | 0.111 |
| 0.10 | Fluor 100 02 20 | 0.130 |
| 0.20 | Fluor 100 02 20 | 0.156 |
| 0.40 | Fluor 100 02 20 | 0.221 |
| 0.59 | Fluor 100 02 20 | 0.279 |
| 0.77 | Fluor 100 02 20 | 0.335 |
| 1.00 | Fluor 100 02 20 | 0.395 |

Calibration Curve

*serial dilution

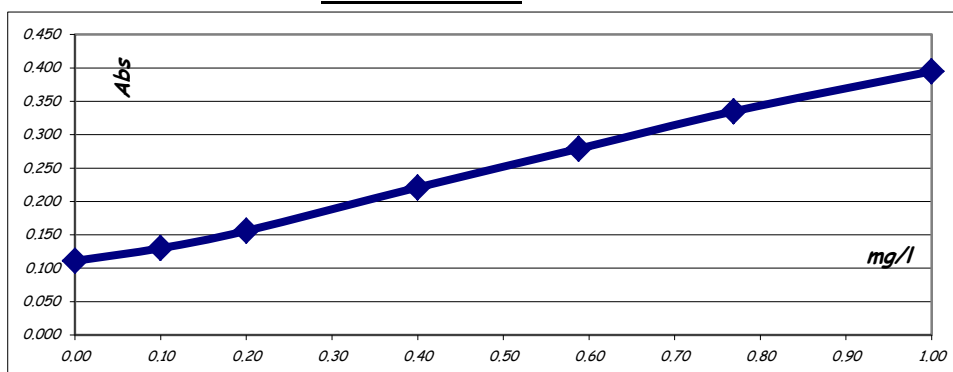

| Test | Result | Specification | Pass / Fail |
|-------------------|---------|----------------|-------------|
| Intercept | 0.111 | 0.1 to 0.116 | Pass |
| Sensitivity | Abs/ppm | | |
| 0 - 0.100 ppm | 0.190 | 0.150 to 0.250 | Pass |
| 0.100 - 0.200 ppm | 0.260 | 0.250 to 0.350 | Pass |
| 0.200 - 0.400 ppm | 0.325 | 0.300 to 0.350 | Pass |
| 0.400 - 0.588 ppm | 0.309 | 0.300 to 0.350 | Pass |
| 0.588 - 0.769 ppm | 0.309 | 0.280 to 0.320 | Pass |
| 0.769 - 1.000 ppm | 0.260 | 0.220 to 0.280 | Pass |

Notes

- Store between 2 & 8 °C
- *Stable on Board Fluoride 1, 2 & 3 5days
- Keep Tightly Sealed when not in use.
- Do not return unused portions to container
- Barcode ID- Fluoride 1 131; Fluoride 2 132; Fluoride 3 133
- * Based on results obtained in our Research Facility

Alan Coonerty
A.Coonerty

An ISO 9001: 2008 Registered Company
Annacotty Susinness Park, Annacotty, Limerick, Ireland
Telephone 061 358190 Facsimile 061 358 191

Hazard Identification

Labelling according to Regulation (EC) No 1272/2008